

PORTABLE PANEPINTO SLING

PORTABLE RESTRAINT



Humane, user-friendly, portable restraint developed in partnership with Panepinto Associates for your micro/minipigs, dogs, goats, lambs and other small to medium sized animals. Exclusively available from Lomir Biomedical, this newest innovation in comfortable restraint is collapsible and lightweight. It can be quickly & easily folded up for trouble-free transport or storage in facilities with limited space or mobile veterinarians. For more information regarding proper use, assembly or hands-on training customized for personnel and facility needs, please visit www.panepinto.com.

1. Place the 4 shorter bars with center welded support wedges marked AB & CD together so the welded mid-sections fit together to form an "X". Make sure the open-hole ends are on top.

2. Connect the AB and CD parts with the long cross bar marked E by inserting the 2 long wing nut bolts into the holes located in the center of each "X" (Figure 1).

3. Tighten the bolts on each end and stand the frame up.

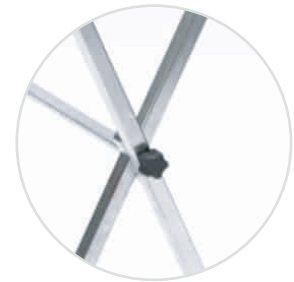


Figure 1



Should you have any questions or need more information please contact our Customer Service Department toll free in North America at 1-877-425-3604 or send an e-mail to info@lomir.com.

Lomir Biomedical Inc.

95 Huot, Notre-Dame-de-l'Île-Perrot, QC
Canada J7W 1M3
Tel: 514-425-3604
Toll Free: 1-877-425-3604
Fax: 514-425-3605
Email: info@lomir.com

Lomir Biomedical Inc.

C/O Agenda Life Sciences
UK PO Box 24 Hull
HU128YK, UK
Tel: 03456 44 55 45
Email: info@lomir.com

Lomir Inc.

PO Box 778
Malone, NY
USA 12953
Tel: 518-483-7697
Toll Free: 1-877-425-3604
Email: info@lomir.com

PORTABLE PANEPINTO SLING

PORTABLE RESTRAINT

- Slide 2 remaining bars (F&G) through sleeves on hammock and place angled ends into holes on top of the "X's" at each end (Figure 2).
- Secure top cross bars on all 4 corners using the clips provided in the shipping box. (Figure 3).
- Place animal into hammock. (Figure 4).
- Use Velcro flaps to secure animal. (Figure 5).
- To remove animal, remove each of the 4 clips which secure the cross bars to the frame and lift the pig while still in the hammock and carry back to pen while contained in the hammock. (This should be done with 2 handlers. If there aren't 2 available, just lift animal out, providing as much body support as possible).
- To fold up for transport or storage, loosen each wing nut and pull the welded joints apart. (Figure 6). Do not remove pin or take completely apart unless planning to disassemble as storage and transport are facilitated if unit stays together.
- Fold up for transport or storage and secure hammock with bungee cord and rope. (Figure 7).



Figure 2



Figure 3



Figure 4



Figure 5



Figure 6



Figure 7

Lomir Biomedical Inc.

95 Huot, Notre-Dame-de-l'Île-Perrot, QC
Canada J7W 1M3
Tel: 514-425-3604
Toll Free: 1-877-425-3604
Fax: 514-425-3605
Email: info@lomir.com

Lomir Biomedical Inc.

C/O Agenda Life Sciences
UK PO Box 24 Hull
HU128YK, UK
Tel: 03456 44 55 45
Email: info@lomir.com

Lomir Inc.

PO Box 778
Malone, NY
USA 12953
Tel: 518-483-7697
Toll Free: 1-877-425-3604
Email: info@lomir.com