

LOMir
BIOMEDICAL INC.

MOBILE RESTRAINT UNIT

Rugged, versatile and humane immobilization of large animals for veterinary and laboratory procedures.

MOBILE RESTRAINT UNIT

A stable immobile and comfortable subject is a key component of any investigation or procedure involving large animals. The innovative Mobile Restraint Unit embodies the most efficient incorporation of all the necessary features. Combined with quality materials and Lomir's attention to detail and workmanship, this device is unique in providing humane immobilization for the animal and efficient, safe, working platforms for the handler. **Special Requirements?** While our inventory of frames and slings will answer the needs of most customers, Lomir is fully equipped to fulfill orders requiring the fabrication of custom-made units for special situations or tasks. For more information, or to order please visit us online at www.lomir.com or call 1-877-425-3604.

FRAME CONSTRUCTION

1. Deep "V" opening in frame for unimpeded access to the head of the animal.
2. Low cross members and deep opening for minimum lifting height and easier positioning of animal.
3. Sling Support bar extenders for lifting. Removable to facilitate access to the animal during procedure and prevent snagging during movement.
4. Available in a range of sizes. Constructed from precision welded 1" square tube stainless steel stock for lightness and strength.



Self-locking, quick release cleats for limb ropes at each corner.



Lower sling support bar brackets for easier access.



Locator pins on all corners to adjust sling support bars.




Corner posts on heavy duty castors with ball-bearing pivots, two with foot activated brakes.



SLING CONSTRUCTION

1. Basic construction of robust, washable canvas type material.
2. Wide reinforced side tubes for easy insertion and removal of adjustable side support bars.
3. All leg and waste apertures edged in cotton jersey for additional comfort.
4. Nylon webbing sewn in at key points throughout for additional strength.

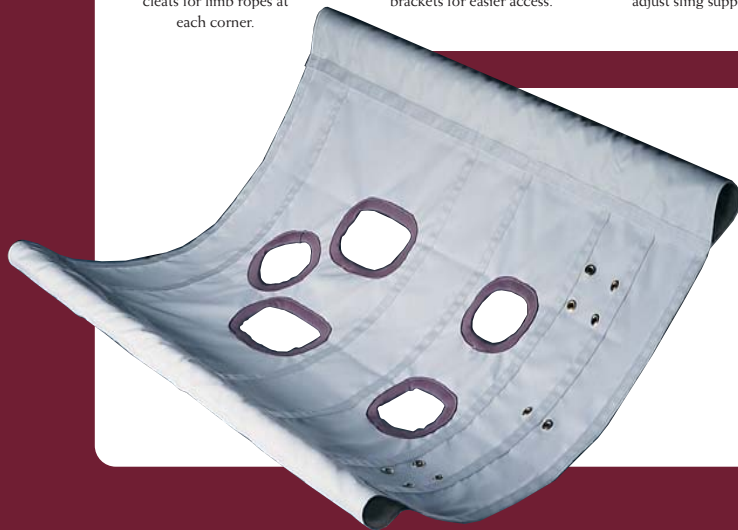

458 East Main Street, Malone, NY, 12953
Tel: (518) 483-7697 Fax: (518) 483-8195
Toll Free: (877) 425-3604
www.lomir.com

95 Huot, Notre-Dame-de-l'Île-Perrot, Qc., J7V 7M4
Tel: (514) 425-3604 Fax: (514) 425-3605
Toll Free: (877) 425-3604
www.lomir.com