



**LOMir**  
BIOMEDICAL INC.



## THE SNUGGLE FAMILY

Protected by U.S. Patent #5,975,028. All Rights Reserved.

Introducing the complete "**Snuggle Family**" the ultimate in stress-free immobilization and versatility for procedures involving guinea pigs, rodents or rabbit subjects.



**"Guinea Pig Snuggle"**



**"Rodent Snuggle"**



**"Bunny Snuggle"**

# THE SNUGGLE FAMILY

Protected by U.S. Patent #5,975,028. All Rights Reserved.

## HANDLE WITH CARE.


Small animals are prone to stress and require **extra care** in handling. Mechanical restraints frequently require the involvement of two personnel: one to restrain and one to perform. The resulting stress and discomfort for the subject can often jeopardize the efficiency of the procedure.

The **"Snuggle Family"** uses the principle of a cocoon rather than a rigid structure for immobilization. The multiple flaps perform two functions: They are adjustable around the subject for effective yet comfortable restraint, and positioned, or left open for access depending on the procedure. Both restraint and procedure can be carried out by one handler.

**Blood Collection, Vaginal Samples, Bolus Injection, Veterinary Exams, Gavage, Transport and more.**

## THE LOMIR DIFFERENCE.

The LOMIR **"Snuggle Family"** is manufactured from flexible yet extremely durable canvas fabric. The fabric is coated with nylon to facilitate laundering for reuse. Wide flaps with velcro fasteners allow for maximum flexibility, ease of use and adjustment for fit between individual subjects. The **"Snuggle Family"** is available in different sizes for a complete range of animals including pregnant females.

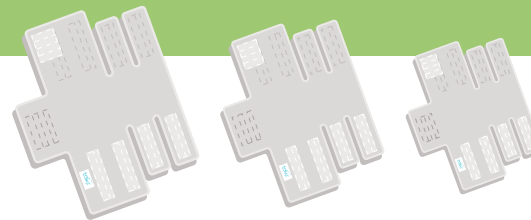

### Sizing Available:

Rabbit - Small, Medium, Large, X-Large

Guinea Pig - Small, Medium, Large, X-Large

Rodent - X-S, Small, Medium, Large, X-Large

458 East Main Street, Malone, NY, 12953  
Tel: (518) 483-7697 Fax: (518) 483-8195  
Toll Free: (877) 425-3604  
www.lomir.com

95 Huot, Notre-Dame-de-l'Île-Perrot, Qc., J7V 7M4  
Tel: (514) 425-3604 Fax: (514) 425-3605  
Toll Free: (877) 425-3604  
www.lomir.com